



```
mne.gui.coregistration(subject=subject,inst=inst)
```

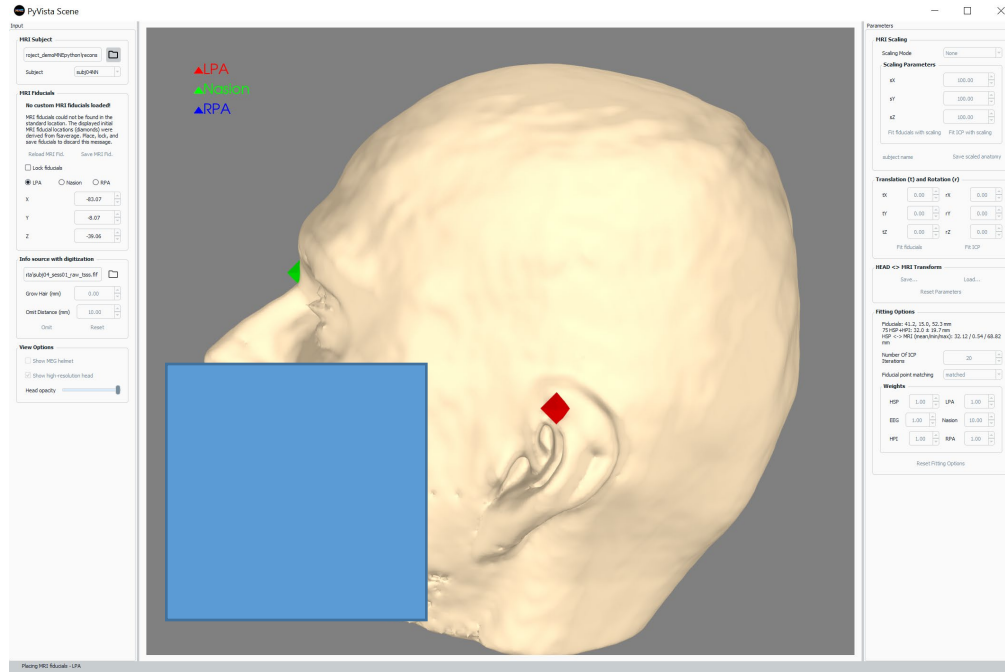

Right click on surface to denote LPA:

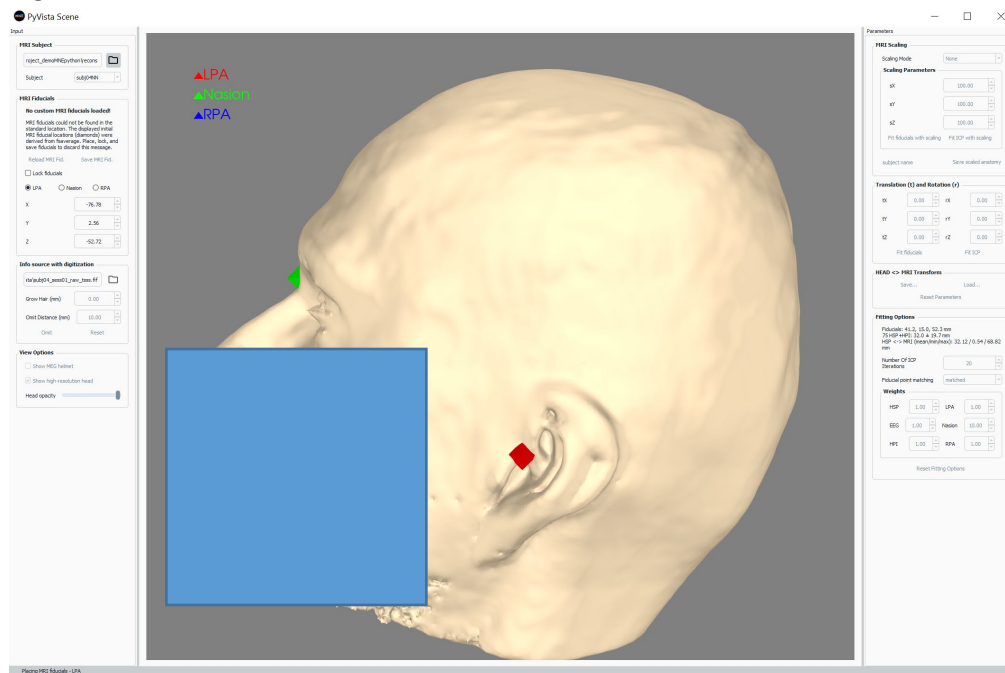

Select 'Nasion' from bullet list:

MRI Fiducials

No custom MRI fiducials loaded!

MRI fiducials could not be found in the standard location. The displayed initial MRI fiducial locations (diamonds) were derived from fsaverage. Place, lock, and save fiducials to discard this message.

Reload MRI Fid. Save MRI Fid.

☐ Lock fiducials


☐ LPA ☒ Nasion ☐ RPA

X-1.55

Y87.29

Z3.37

Right click on surface to denote Nasion:

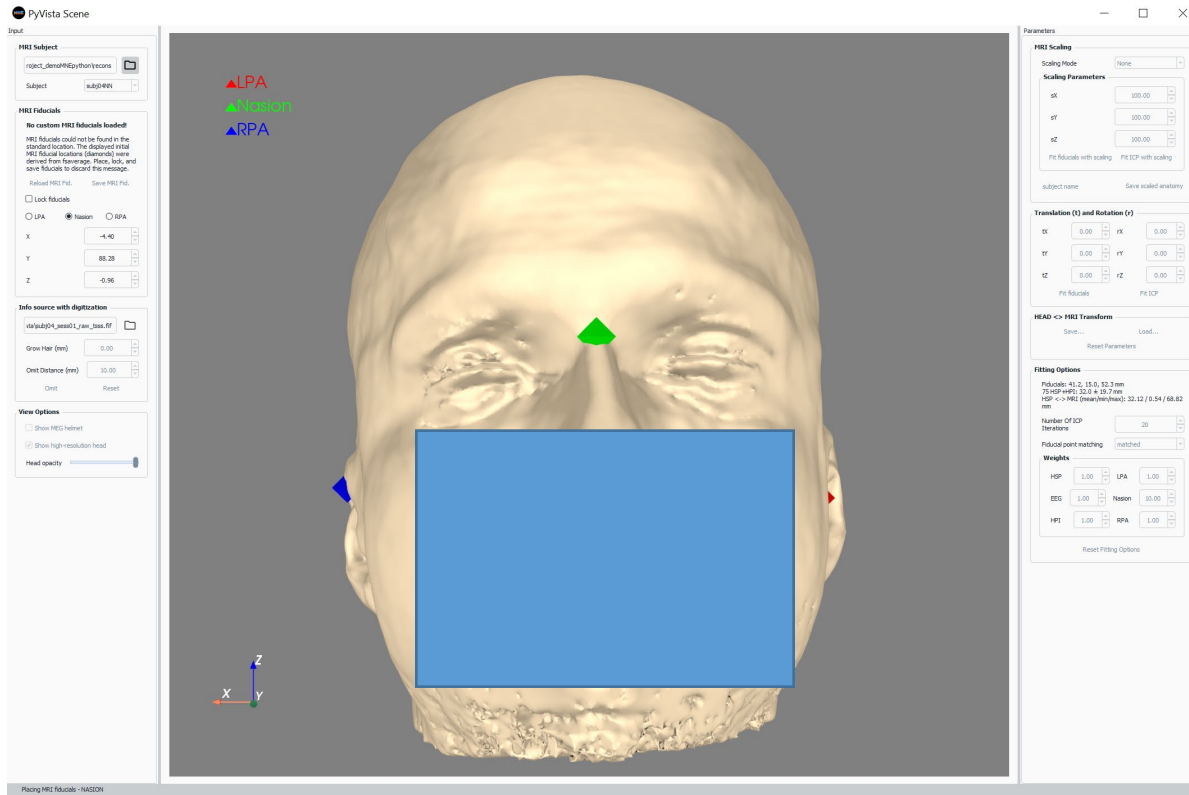

Select 'RPA from bullet list:

MRI Fiducials

No custom MRI fiducials loaded!

MRI fiducials could not be found in the standard location. The displayed initial MRI fiducial locations (diamonds) were derived from fsaverage. Place, lock, and save fiducials to discard this message.

Reload MRI Fid.

Save MRI Fid.

☐ Lock fiducials

☐ LPA

☐ Nasion

☒ RPA

X

75.88


Y

-8.72

Z

-50.56

Right click on surface to denote RPA:

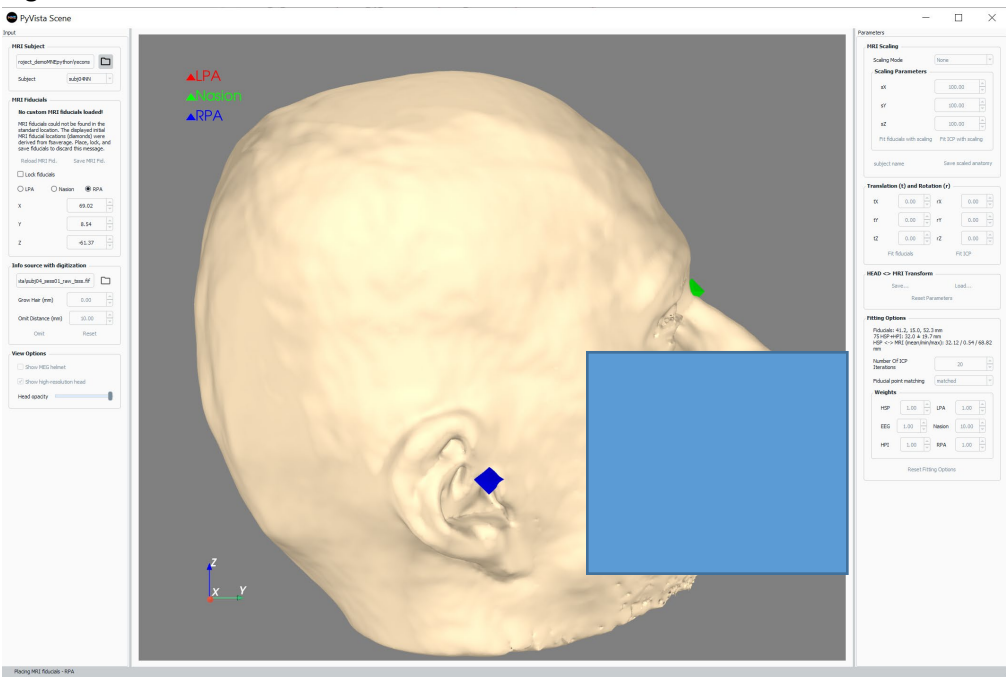

Lock fiducials:

MRI Fiducials

No custom MRI fiducials loaded!

MRI fiducials could not be found in the standard location. The displayed initial MRI fiducial locations (diamonds) were derived from fsaverage. Place, lock, and save fiducials to discard this message.

Reload MRI Fid. **Save MRI Fid.**

☒ Lock fiducials

☐ LPA ☐ Nasion ☒ RPA

X 69.02

Y 8.54

Z -61.37

Click 'Save MRI Fid.' This will create a <subject>-fiducials.fif file inside the bem folder of the subject freesurfer directory.

MRI Fiducials

No custom MRI fiducials loaded!

MRI fiducials could not be found in the standard location. The displayed initial MRI fiducial locations (diamonds) were derived from fsaverage. Place, lock, and save fiducials to discard this message.

Reload MRI Fid. **Save MRI Fid.**

☒ Lock fiducials

☐ LPA ☐ Nasion ☒ RPA

X 69.02

Y 8.54

Z -61.37

Press Fit fiducials. Then press Fit ICP (a few times, as appropriate to refine registration)

Translation (t) and Rotation (r)

tX 0.00 rX 0.00

tY 0.00 rY 0.00

tZ 0.00 rZ 0.00

Fit fiducials **Fit ICP**

